

Revised 10 mg Pharmacist Card

Diastat AcuDial™ (diazepam rectal gel) IV
10 mg Delivery System


0 01870 65820 9
0187-0658-20

ATTENTION PHARMACIST

- ADJUST SYRINGES PRIOR TO DISPENSING**
- Follow instructions on back of this card to set prescribed dose



PHARMACIST MUST
DIAL AND LOCK
CORRECT PRESCRIBED DOSE

- Product should not be dispensed until **BOTH** syringes have been dialed and locked
- INSTRUCT CAREGIVERS/PATIENTS** to confirm that the correct dose has been dialed and review instructions on bottom of tray before they leave the pharmacy
- Do not** place prescription label on front panel until **this card has been removed**


128763

Available Doses from this delivery system

5 mg7.5 mg10 mg

Contains two pre-filled, **ADJUSTABLE** dose syringes for rectal delivery with lubricating jelly and instructions for use


DIASTAT® AcuDial™ PHARMACIST INSTRUCTIONS

— READ BEFORE COMPLETING —

1. HOLD

Remove a syringe from case, hold barrel of syringe as shown.

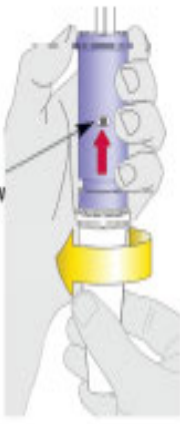
Do Not Remove Cap!



Labels: Barrel, Cap

2. ADJUST


Grasp cap firmly with other hand and turn to adjust dose.



Labels: Dose Window

3. LOCK

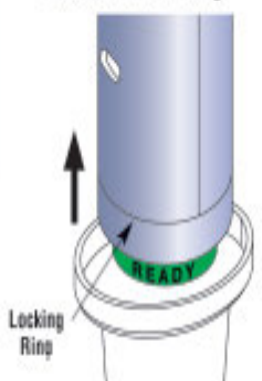CONFIRM prescribed dose shows in window. Now grasp locking ring and push upward to lock both sides of ring.



Labels: Locking Ring


4. REPEAT

Repeat steps 1-3 for second syringe. Return both syringes to case.


DISCARD
THIS CARD
WHEN
FINISHED

FOR QUESTIONS, OR IF RING IS LOCKED AT INCORRECT DOSE, CALL 877-361-2719.

Revised 20 mg Pharmacist Card

Diastat AcuDial 
(diazepam rectal gel)
20 mg Delivery System


VALEANT
Pharmaceutical North America



0 01870 65920 6
0187-0659-20

ATTENTION PHARMACIST

- ADJUST SYRINGES PRIOR TO DISPENSING**
- Follow instructions on back of this card to set prescribed dose



PHARMACIST MUST
DIAL AND LOCK
CORRECT PRESCRIBED DOSE

- Product should not be dispensed until **BOTH** syringes have been dialed and locked
- INSTRUCT CAREGIVERS/PATIENTS** to confirm that the **correct dose** has been dialed and review instructions on bottom of tray **before** they leave the pharmacy
- Do not** place prescription label on front panel until this card has been removed



128764

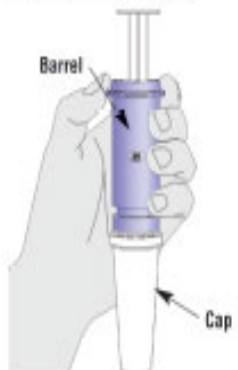

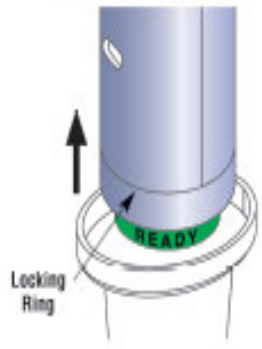
Available Doses from this delivery system

12.5 mg	15 mg	17.5 mg	20 mg
---------	-------	---------	-------

Contains two pre-filled, ADJUSTABLE dose syringes for rectal delivery with lubricating jelly and instructions for use

DIASTAT® AcuDial™ PHARMACIST INSTRUCTIONS

— READ BEFORE COMPLETING —

- 1. HOLD**
Remove a syringe from case, hold barrel of syringe as shown.
Do Not Remove Cap!

- 2. ADJUST**
Grasp cap firmly with other hand and turn to adjust dose.

- 3. LOCK**
CONFIRM prescribed dose shows in window.
Now grasp locking ring and push upward to lock both sides of ring.

- 4. REPEAT**
Repeat steps 1-3 for second syringe.
Return both syringes to case.

DISCARD THIS CARD WHEN FINISHED

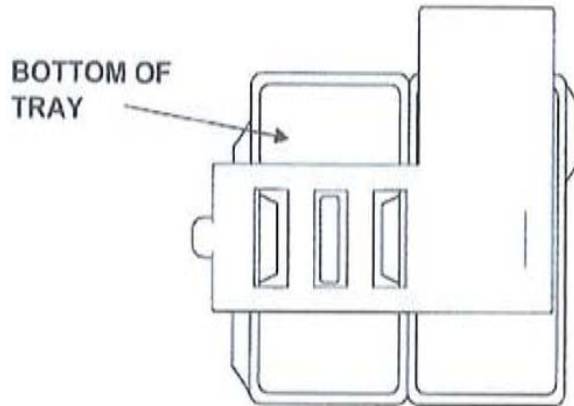
FOR QUESTIONS, OR IF RING IS LOCKED AT INCORRECT DOSE, CALL 877-361-2719.

Assembly Instruction Diagram (non-adhesive band/tab label)

Attachment

ASSEMBLY INSTRUCTIONS:

41. Place one (1) flap insert print side down into thermoformed tray as shown:



42. Place two (2) shroud assemblies into thermoformed tray as shown below, on top of the flap insert. Rotate shroud assemblies such that Lot Code and Expiration Date of each shroud assembly are visible through the side windows at the bottom of the tray. Wrap vertical flap around shroud assemblies and insert tab through the slot in the flap as shown.

